

Century™ Max T5200

Bariatrics Long Term Care Bed

COLOR: ANTIQUE WHITE GLOSS FRAME WITH
LIGHT OAK HEAD & FOOT BOARDS



The Century™ Max long term care bed is a true bariatrics bed which will accommodate a weight limit up to 1,000 lbs. It is specially engineered and designed for large residents or others with chronic impairments that reduce their degree of independence in their daily life's activities. Top manufacturing technology with the highest quality components make these beds the best choice for bariatric users. It's fully electric and the 5 function hand pendent control performs all functions of the bed. It comes in frame lengths of 80" or 84" and bed widths of 48", 54" or 60".

Tuffcare®
HIGHER STANDARD FOR HEALTH CARE™

Century™ Max T5200

Bariatrics Long Term Care Bed

**COLOR: ANTIQUE WHITE GLOSS FRAME WITH
LIGHT OAK HEAD & FOOT BOARDS**



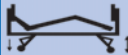


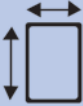



STANDARD FEATURES

- * Welded angle steel construction with a low center reinforced bridge provides maximum stability.
- * Aluminum alloy bed surface has the same strength as steel but weighs 30% less.
- * Two piece split frame design is easy to transport, lighter to set up, and operate.
- * Trendelenburg position, and its reverse, can be easily achieved with the hand control pendant.
- * Built in back up battery can fully operate the bed for 2 hours in the case of a power failure.
- * LED night light is a great for checking on the patient as well as safely guiding them to their bed.
- * Four 5" locking swivel casters provide increased security.
- * Accommodates user weight up to 1,000 lb
- * HCPCS Code: E0304

OPTIONAL ACCESSORIES

- * One pair or two pairs of half rails
- * High density foam mattress
- * Alternating air pad system
- * 4" bed length extender
- * Bed transporter & storage unit
- * Trapeze bar
- * Gel overlay mattress
- * Deluxe convoluted top foam mattress
- * Alternating & low air loss mattress system
- * Deluxe overbed table
- * Deluxe bedside cabinet
- * Patient Lift

SPECIFICATIONS

MODEL	HI/LO OPERATION	HEAD/FOOT OPERATION	FABRIC HEIGHT	OVERALL WIDTH	OVERALL LENGTH	MATTRESS SIZE	HEAD FRAME WEIGHT	FOOT FRAME WEIGHT	BED ENDS WEIGHT
									
T5248	ELECTRIC	ELECTRIC	15" to 28"	48"	86"	48" X 80"	139 lbs.	129lbs	83lbs
T5254	ELECTRIC	ELECTRIC	15" to 28"	54"	86"	54" X 80"	143 lbs.	135lbs	86lbs
T5260	ELECTRIC	ELECTRIC	15" to 28"	60"	86"	60" X 80"	145 lbs.	137lbs	88lbs



Use the rising head and foot surface functions together with the Reverse Trendelenburg function to achieve a lounge chair position.



Offset Leg design ensures the 4 casters stay stationary and that the bed surface does not shift while elevating or descending.



Five function ergonomic hand pendant with LED night light button.

Tuffcare®
HIGHER STANDARD FOR HEALTH CARE™